



HUAYI ISOTOPES

HIC Quality Products are proudly distributed by:

ISO FLEX USA

P.O. Box 29475
San Francisco, CA 94129
415-440-4433
www.isoflex.com

Huayi Isotopes (China)

Phone : +86-512-5249-1118 E-mail : sales@huayi-isotopes.com

CONTENT

I. PET

1. PET DISPENSING SYSTEM	(1)
2. PET Dispensing Pig 30ml	(1)
3. PET Dispensing Pig 10ml	(2)
4. PET Dispensing System Syringe Shield 3/5cc	(2)
5. PET Dispensing System Syringe Shield 10cc	(2)
6. PETPig™ Syringe Pig/Syringe Shield 3/5cc,10cc	(3)
7. PETPig™* Cradle	(3)
8. PET Hot Cell Pig	(3)
9. PET QC Syringe carrier 1cc	(4)
10. Nordion 3ml V-Vial Adapter.	(4)
PET Vial Adapter 10ml	(4)
PET Vial Adapter 10ml-DuPont	(4)
PET Vial Adapter 20ml	(4)
PET Vial Adapter 25ml	(4)
Vial Adapter For PET Hot Cell Pig	(4)
11. PET L-Block Shield	(4)
12. PET Clean Air Hood	(5)
13. PET Waste Cart	(6)
14. PET Charcoal Filter (Lg:)	(6)
15. Waste Box for gloves	(6)
16. PET Body Shield	(6)
17. PDP Transportation Cart	(6)
18. PET QC Vial Holder.....	(6)
19. FDGPig for 3/5cc Syringe.....	(6)
20. PET Dispensing System Syringe Shield 3/5cc With Lead Glass Window.....	(7)
21. PET Lead Trans Pig.....	(7)

II. Tech

1. Single Product Dispensing Station	(8)
Multi-Product Dispensing Station	(8)
Tungsten Magnetic Vial Shield	(8)
2. The EluTer Elution Vial Shield for DuPont Technelite Generator	(9)
3. SECURE+™ Unit Dose Syringe Pig	(9)
SECURE+™ Unit Dose Syringe Pig Tray	(9)
SECURE+™ Unit Dose Syringe Pig 10cc	(9)
4. SECURE® Bio Cap Dispenser	(10)
5. Hamilton Vial Heating Shield	(10)
6. Defender™ Dose Drawing Syringe Shields	(11)
7. Lifetime™ Syringe Shields	(12)
Latchkey Syringe Shield™	(12)
8. SECURE® Syringe Shield	(13)
9. Latchkey Beta/Gamma Syringe Shield 10cc	(13)
10. Centrifuge Tube Shield Cup	(13)
11. Tech L-Block.	(14)
12. QC Station	(14)
13. QT L-Block	(14)
14. Tech Clean Air Hood	(15)
15. Tech Lead Waste Barrel (0.75cm)	(16)
Tech Lead Waste Barrel (1.27cm)	(16)
16. Generator Box	(16)
17. Clean Air Hood (Small Type).....	(17)
18. The EluTer Elution Vial Shield for DuPont Technelite Generator(Type 2).....	(17)
19. The EluTer Elution Vial Shield for DuPont Technelite Generator(Type 3).....	(18)
20. R-Lead Unit Dose Syringe Pig.....	(18)
21. R-Lead Magnetic Vial Shield.....	(18)
22. Safe Tview Syringe Shield 1cc-BD 1ml Luer Luck.....	(19)
Safe Tview Syringe Shield 2cc-Braun 2ml	(19)
Lifetime Syringe Shield 1cc-BD 1ml Luer Luck.....	(19)
Lifetime Syringe Shield 1cc- Braun 2ml.....	(19)
23. R-Lead Centrifuge Tube Shield Cup.....	(19)

III. Iodine

1. Iodine Dispensing System.....	(20)
2. Piglet™.....	(20)
Piglet™ Sling.....	(20)
3. Piglet 2™.....	(21)
Piglet™ Pooper.....	(21)
Gotti Cup.....	(21)
4. Piglet 3™.....	(21)
5. SECUREplus™ Delivery Multi Dose Nylon Case Foam Insert-one Hole.....	(21)
SECUREplus™ Delivery Multi Dose Nylon Case Foam Insert-two Holes.....	(22)
6. Seed Dipper.....	(22)
7. Iodine Charcoal Filter Box (Sm.).....	(22)
Charcoal Air Filer Type 787 (I-131) For Iodine Charcoal Filter Box & Glove Box.....	(22)
8. I-131 Waste Storage Cabinet.....	(22)
9. Iodine Lead Waste Barrel (1.9cm thick).....	(23)
10. I-131-L -Block.....	(24)
11. Glove Box & Table.....	(24)
12. NeRD.....	(25)
13. Tweezers - 9.75 Inch.....	(25)
14. Forceps 10", Straight.....	(25)
Forceps 10", Curved.....	(25)
Forceps Guard(round).....	(25)
Forceps Guard (trapezoid).....	(25)
15. R-Iodine Pig 1.....	(26)
16. R-Iodine Pig 2.....	(26)
17. Squeezer QC Syringe Shield.....	(26)
18. Lead Loose Seed Holder.....	(26)

IV Stainless Steel

1. SS Cabinet:(A)-24" Drawers.....	(27)
SS Cabinet:(C)-36" Cabinet.....	(27)
SS Cabinet:(C)-48" Cabinet.....	(27)
SS Cabinet:(D)-48" Sink.....	(27)
2. SS Cabinet:(E)-24" Overhead Cabinet.....	(27)
SS Cabinet:(F)-48" Overhead Cabinet.....	(27)
3. SS Cabinet:(G)-48" Overhead Open Shelf.....	(28)
4. Conveyer Table.....	(28)
5. Table Base for 2 Waste Barrels.....	(28)
Table Base for 2 Waste Barrels(Special).....	(28)
6. Table Base for 3 Waste Barrels.....	(29)
7. Lead Funnel & Tube.....	(29)
8. SS Cabinet:(SL)-Right Support Leg.....	(29)
SS Cabinet:(SL)-Left Support Leg.....	(29)
9. SS Cabinet:(CG)-Corner Guard.....	(29)


V. Others

1. Ashley Dose Calibrator Funnel Dipper.....	(30)
Standard Dose Calibrator Dipper.....	(30)
Well Liner Sleeve.....	(30)
2. Hot Box.....	(31)
3. SS WBC Case.....	(31)
4. Seed Mick Cartridge Holder.....	(32)
Loose Seed Holder.....	(32)
5. Sing Core Filter.....	(33)
6. J-Tungsten Syringe Shield.....	(33)
J-Tungsten Stopper Shield For Cup Shield.....	(33)
7. Tungsten Funnel.....	(33)
8. W-Cu Alloy Hollow Tube.....	(33)
W-Cu Alloy Solid Bar.....	(33)

PET

I. PET DISPENSING SYSTEM*

The solid tungsten PET Dispensing System* permits the safe dispensing of high-energy PET radiopharmaceuticals without the need for expensive remote handling systems. The magnetic "docking" feature of the solid tungsten syringe shield into the dispensing pig allows low exposure dose drawing. It is designed with 3.5 cm thick walls and accommodates up to a 30ml vial.

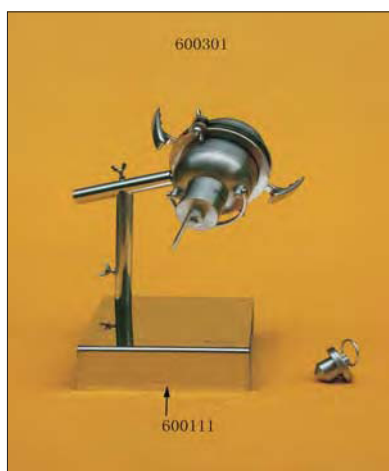

The 600301 PET Dispensing System* includes:

- 600111 Stainless Steel Base , Armature & Ring
(抽药用不锈钢支架)
316×210×60mm N.W:5Kg
- 600304 PET Dispensing Pig*30ml
(PET药液存放罐 30ml)
Ø108.5×157.5mm N.W:20.5Kg
- 600319 PET Dispensing Pig10ml
(PET药液存放罐 10ml)
Ø97.5×144mm N.W:16Kg
- 600302 PET Dispensing Syringe Shield 3cc & 5cc*
(PET3cc/5cc针管防护套)
Ø57×85mm N.W:2.6Kg
- 600303 PET Dispensing System Syringe Shield 10cc
(PET10cc针管防护套)
Ø59×98.5mm N.W:3.4Kg

II. PET Dispensing Pig 30ml

The PET Dispensing Pig provides a full 3.5 cm of solid tungsten to safely shield high energy PET radiopharmaceuticals . Designed to accommodate 30ml vials (adapters are available for smaller vials). The threaded lid attaches to the body of the pig so that no radiation "shine through" occurs. The tungsten stopper, with attached pull ring ,is held in place magnetically. The handle permits the pig to be easily carried as a transport container as well as a dispensing pig. Weight 43 lbs.

- 600304 PET Dispensing Pig*30ml
(PET药液存放罐 30ml)
Ø108×157.5mm N.W:20.5Kg
- 600311 PD Pig Ammo Can Foam Insert 30ml
(药液存放罐上下海绵)
280×140×115mm N.W:0.1Kg
- 600312 PDS Circular Sponge (100/bag)
(PET药液存放罐内黄海绵片)
Ø35.8mm
N.W:0.04Kg/50pcs
- 600313 PET Dispensing Ring
(PET30ml药液存放罐套环)
Ø138×Ø109×5mm N.W:0.82Kg
- 600326 PET Vial Adapter 10ml
(PET30ml药液存放罐内衬杯)
Ø36×62mm N.W:0.05Kg



III. PET Dispensing Pig 10ml

- 600319 PET Dispensing Pig 10ml
(PET 10ml 药液存放罐)
Ø 97.5 × 144mm N.W: 16Kg
- 600320 PET Armature & Ring 10cc
(PET 10ml 药液存放罐网拍)
Ø 127 × Ø 98 × 5mm N.W: 0.77Kg
- 600322 PDS Ammo Can Foam Insert 10cc
(PET 10ml 药液存放罐海绵)
280 × 140 × 115mm N.W: 0.1Kg



IV. PET Dispensing System Syringe Shield 3/5cc

The PET Syringe Shield magnetically docks with the PET Dispensing Pig. Designed to accept 3cc and 5cc B-D syringes, it places the needle inside the vial septum when engaged. The body of the shield is constructed of solid tungsten 2.1cm thick for maximum shielding. The external calibration rod allows the precise volume to be withdrawn without a leaded glass viewing port, where high exposure levels cannot be adequately shielded. The system allows the plunger to be pulled back with a pair of forceps, allowing you to keep your exposure ALARA. Weight 5.8 lbs

- 600302 PET Dispensing System Syringe Shield 3/5cc
(PET 3cc/5cc 针管防护套)
Ø 57 × 85mm N.W: 2.6Kg
- 600316 PET: 3cc Gauge
(PET 3cc 针管防护套剂量标尺)
Ø 23 × 88.5mm N.W: 0.06Kg
- 600317 PET: 5cc Gauge
(PET 5cc 针管防护套剂量标尺)
Ø 24 × 78.5mm N.W: 0.05Kg
- 600332 PET Dispensing Syringe Shield Holder (short)
(PET 3cc-5cc 针头防护套)
Ø 64 × 55mm N.W: 1.77Kg

V. PET Dispensing System Syringe Shield 10cc

- 600303 PET Dispensing System Syringe Shield 10cc
(PET 10cc 针管防护套)
Ø 59 × 98.5mm N.W: 3.4Kg
- 600318 PET: 10cc Gauge
(PET 10cc 针管防护套剂量标尺)
Ø 26 × 100.5mm N.W: 0.07Kg
- 600339 PET Dispensing Syringe shield Holder (long)
(PET 10cc 针头防护套)
Ø 64 × 75mm N.W: 2.24Kg



VI. PETPig™ Syringe Pig/Syringe Shield

PETPig™ permits the safe transport and administration of unit dose PET radiopharmaceuticals. It's constructed with solid tungsten walls 2.25cm thick and yet weighs only 15.6lbs. The "T" handle on the PETPig™ cap allows the unit to be easily lifted out of traditional "ammo can" delivery cases. The threaded top provides easy access to the 3 or 5cc syringe. The use of the thermos style handle reduces hand exposure by permitting the PETPig™ to be carried to the imaging suite without holding container sidewalls. Prior to injection, the base unscrews, allowing the center portion to be used as a syringe shield. When placed in the optional PETPig™ Cradle (Product #600308), patient administrations can be performed with ease and maximum shielding.



- 600300 PETPig™ Syringe Pig/Syringe Shield 3/5cc*
(PET3cc/5cc药物针管运输罐)
Ø 60×241mm N.W: 7.5Kg
- 600309 PETPig™ Syringe Pig/Syringe Shield 10cc*
(PET10cc药物针管运输罐)
Ø 63.5×244.7mm N.W: 9.5Kg
- 600333 PETPig 3/5cc Sleeve
(PET3cc/5cc药物针管运输罐底座防护套)
Ø 78×Ø 66×123mm N.W: 1.74Kg
- 600307 PETPig™ Ammo Can Foam Insert 3/5cc
(PET3cc/5cc药物针管运输罐铁箱用海绵)
280×140×120mm N.W: 0.08Kg
- 600321 PETPig™ Ammo Can Foam Insert 10cc
(PET10cc药物针管运输罐铁箱用海绵)
280×140×120mm N.W: 0.11Kg
- 600310 PETPig™ Donut Sponge (100/bag)
(PET3cc/5cc药物针管运输罐黄海绵环)
Ø 35.8×Ø 15mm
N.W: 0.04Kg/50pcs

VII. PETPig™ Cradle

The PETPig™ Cradle holds the PETPig™ 3cc and 5cc syringe pig/syringe shield* in place so you can unscrew the top and bottom portions. It holds the shielded syringe securely and permits easy, quick injections.

- 600308 PETPig™ Cradle 3/5cc
(PET3cc/5cc药物针管运输罐不锈钢支架)
89.5×112×68mm
N.W: 0.6Kg
- 600325 PETPig™ Cradle 10cc
(PET10cc药物针管运输罐不锈钢支架)
110×115×68mm
N.W: 0.85Kg

VIII. PET Hot Cell Pig

- 600329 PET Hot Cell Pig
(PET高剂量药物存放罐)
Ø 136×135mm N.W: 39Kg
- 680003 Vial Adapter For PET Hot Cell Pig
(PET高剂量药物存放罐内衬杯)
Ø 56×87mm N.W: 0.11Kg




IX.

- 600347 PET QC Syringe carrier 1cc**
(PET1cc质控存放罐)
Ø52×108.5mm N.W:3.7Kg



X.

- 600132 Nordion 3ml V-Vial Adapter (I¹³¹ 药物存放罐内衬杯)**
Ø25.5×48mm N.W:0.02Kg
- 600326 PET Vial Adapter 10ml (PET药液存放罐内衬杯10ml)**
Ø36×62mm N.W:0.05Kg
- 600327 PET Vial Adapter 10ml-DuPont (杜邦PET药液存放罐内衬杯)**
Ø35.5×62mm N.W:0.04Kg
- 600328 PET Vial Adapter 20ml (PET药液存放罐内衬杯20ml)**
Ø36×62mm N.W:0.04Kg
- 600336 PET Vial Adapter 25ml (PET药液存放罐内衬杯25ml)**
Ø35.5×62mm N.W:0.02Kg
- 680003 Vial Adapter For PET Hot Cell Pig (PET高剂量药物存放罐内衬杯)**
Ø56×87mm N.W:0.11Kg

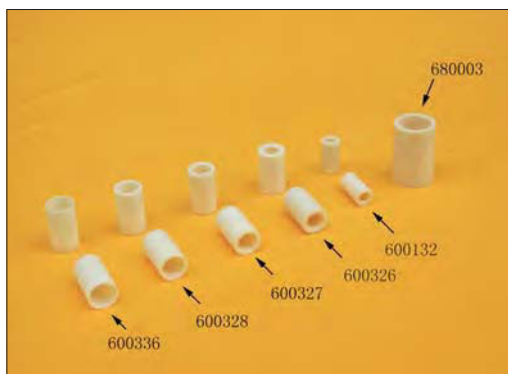

XI. PET L-Block

Specially designed for PET facilities or nuclear medicine departments that are working with high-energy positron-emitting isotopes. Sturdy steel PET table makes rolling and maneuvering exceptionally easy. The ideal addition to labs performing these specialized exams.

- 600314 PET L-Block Shield**
(PET药物分装 L型防护屏)
Dimensions: 380×365×600mm
Lead: thick 65mm
Lead Glass: 225×145×120mm
Material: Lead, Stainless steel
N.W: 270Kg

- 600315 PET Stainless Steel Table**
(PET药物分装 L型防护屏移动车)
Dimensions: 915×700×820mm
N.W: 52Kg



XII.PET Clean Air Hood

This drawing station made of easy-to-clean stainless steel, features a built-in well for the dose calibrator, custom-sized shielded wells for waste collection, a fully shielded base and an integral L-Block style shield. In addition, the drawing station features a 0.3 micron filter, built-in lighting and a pulley and cable assembly system to raise and lower syringes from the ion chamber. This device is the perfect solution for those facilities that use large amounts of activity and prepare doses for many studies.

SPECIFICATIONS:

Working Area: 1100×710×740mm Lead Shielded

Color: Stainless Steel

Finish: Smooth

Electrical: 115V, 60Hz, single phase, 15 amps, 9' power cord, 4 wire

Final Filter: 18.00×48.00×3.00 Aluminum Frame, 99.909% efficient, .3 microns w/neoprene air seal gasket w/spring tension, spring tension compression clamps.

Total CFM: 600

Lamps: 3-F48T12CW fluorescent cool white, slimline.

Ballast: F36/48T Slimline, 115V, 60Hz, single phase

Motor: 0.5"hp, 1200 rpm

Control: Solid state electronic motor speed control

Material: 16 gauge cold rolled steel

Work Top: 0.5" formed stainless steel front/rear return panels

Outlet: 115V, 60Hz, single phase, 5 amps, grounded, universal jack

Hoist: Pulley and cable assembly for inserting and removing product from ion chamber

Waste Chute: Lead shielded 8" I.D × 12" deep cylinder

Shield: Stainless steel personnel shield lead lined bottom and picture frame with 1" lead glass

Sleeve: Stainless steel formed calibrator holder

Slide: Stainless steel calibrator slide rail, adjustable for viewing

TABLE STAND

Size: 1280×840×980mm

Color: Stainless Steel

Finish: Smooth

Material: 11 gauge, 2×2 cold rolled steel tube and 16 gauge 1×1 cold rolled steel

Stand with Wheels: moveable or stand

Size: 1280×840×2540mm

Lead: thick 12mm

Lead Galss: 510×380×40mm

Material: stainless steel thick 2mm

Weight: 720kg



600330

PET Clean Air Hood

(PET药物分装100级无菌净化通风柜)

XIII. PET Waste Cart

The perfect accessory for departments using PET radionuclides. This cart incorporates the convenience of a fully shielded well that is on a movable cart. The well is located within a lockable storage compartment that serves to further isolate the waste.

600334 PET Waste Cart
(PET废弃物存放移动车)
950×485×810mm N.W: 77Kg
Material: Lead, Stainless Steel



XVI.



600338 PET Body Shield
(PET药物制取设备
维修移动防护屏)
Size:
650×550×1727mm
Lead: thick 40mm
Lead Glass:
240×165×120mm
Material:
Lead, Stainless Steel
N.W: 449 Kg

XVII.

600335 PDP Transportation Cart
(PET药物存放罐
移动车)
Ø610×860mm
Material: Stainless Steel
N.W: 13.5Kg



XIV. 600324

PET Charcoal Filter (Lg.)
(尘埃吸附式过滤过滤器)
Size: 950×680×715mm
Material: Stainless Steel
N.W: 161Kg

600345

Filter Core
(PET放射性药物制取室
双芯过滤器滤芯)
Size: 600×290×600mm
N.W: 45Kg



XV. 600341

Waste Box for gloves
(PET药物作废手套箱)
365×265×827mm
Lead: 15mm
Material: Lead
Stainless Steel
N.W: 45Kg



XVIII.

600363 PET QC Vial Holder
(PET质控针管防护灌)
Ø98.5×100mm
Material: Tungsten
N.W: 11.9Kg



XIX.FDGPig for 3/5cc Syringe

The FDG Pig and the FDG Pig Holders are made of solid tungsten. The FDG Pig is used for putting the syringe into it after dose drawing. And it can be easily lifted out of the FDG Pig Holder without hand exposure.

The shipping bag is made of Nylon. There are two kinds of shipping bag. One has built-in roller wheels and a retractable handle making it very easy to move around. Another is also available without roller wheels and a retractable handle upon request. Included in the shipping bag is a protective, high-density foam compassing the FDG Pig Holder on all sides.



690300 FDG Pig for 3/5 cc. Syringe
(R-PET3/5cc药物运输罐)

Material: Tungsten Size: $\varnothing 44 \times 212\text{mm}$
N.W: 3.5Kg



690302



690302-2

690301 FDG Pig Holder for 3/5 cc. FDG Pig
(R-PET3/5cc药物运输罐护套)

Material: Tungsten Size: $\varnothing 80 \times 228\text{mm}$
N.W: 9.23Kg

690303 Protective Foam for one FDG Pig (hard/soft)
(R-尼龙包用一孔海绵)(硬/软)

Size: $260 \times 160 \times 345\text{mm}$

690304 Protective Foam for two FDG Pig (hard/soft)
(R-尼龙包用二孔海绵)(硬/软)

Size: $260 \times 160 \times 345\text{mm}$

690302 FDG Shipping Bag with rolling wheels and pull-handle
(R-拉杆式尼龙包)

Material: Nylon and PVC

690302-2 FDG Shipping Bag without rolling wheels and pull-handle
(R-手提式尼龙包)

Material: Nylon and PVC



690300 690301



690300 690301

XX.

600323 PET Dispensing System Syringe Shield
3/5cc With Lead Glass Window

(PET3/5cc带可视窗圆桶型防护套)

Material: Tungsten Size: $\varnothing 31 \times 71\text{mm}$
N.W: 0.52Kg



XXI.

600306 PET Lead Trans Pig
(PET10ml铅运输罐)

Material: Lead
Size: $\varnothing 103 \times 135\text{mm}$
N.W: 11.4Kg



600306

Tech

I.

1. Single Product Dispensing Station

The Single Product Dispensing Station permits dispensing of radiopharmaceuticals without having to hold the vial shield. Adjustable stainless steel base provides flexibility and stability. Fully adjustable armature permits positioning to any height or angle. Accepts Tungsten Vial Shield(#600205).Vial shield not included.

2. Multi-Product Dispensing Station

This rotating carousel is designed to hold 6 tungsten vial shields with magnetic lids, permitting the easy withdrawal of patient doses. Once the morning technologist has made all of the required kits, anyone can select the proper produce for each patient scan by simply spinning the carousel. Stainless Steel base accommodates standard lead bricks. The vertical stainless steel pole has multi-level height adjustments. Accepts Tungsten Vial Shield or Tungsten Vial Shield with Swivel Lid vial shields (not included).

3. Tungsten Magnetic Vial Shield

This 6 mm thick tungsten vial shield assures safe shielding of vials containing radioactive liquids. The shield has a removable screw top with a magnetic cap for easy containment of vials. The strong tungsten will keep its shape under the toughest handling is essentially indestructible. Fits most 10 ml vials.



- 600221 **Cardiolite® BoosterSeat**
(⁹⁹Tc药液短瓶存取罐内有机玻璃垫)
Ø 25×17mm N.W:0.01Kg
- 600267 **Tungsten Tech Coaster**
(⁹⁹Tc药液短瓶存取罐周转盘)
Ø 70×10mm N.W:0.57Kg

- 600205 **Tungsten Vial Shield long**
(⁹⁹Tc药液长瓶存取罐)
Ø 43×76mm N.W:1.41Kg
- 600208 **Tungsten Magnetic Vial Shield**
(⁹⁹Tc药液短瓶存取罐)
Ø 39×56.5mm N.W:0.95Kg
- 600269 **Tungsten Magnetic Vial Shield Top**
(⁹⁹Tc药液瓶存取罐上盖)
Ø 39×7mm N.W:0.12Kg
- 600270 **Tungsten Magnetic Vial Shield Middle**
(⁹⁹Tc药液瓶存取罐中盖)
Ø 39×14.5mm N.W:0.14Kg
- 600209 **Single Product Dispensing Station**
(⁹⁹Tc药物抽药不锈钢单头支架)
215×110mm Stainless Steel
N.W:1.51Kg
- 600210 **Multi-Product Dispensing Station**
(⁹⁹Tc药物抽药不锈钢多头支架)
215×110mm Stainless Steel

II. The EluTer Elution Vial Shield for DuPont Technelite® Generator

Solid tungsten vial shield designed to shield multi-curie elutions. The tough tungsten will retain its shape under the roughest handling and is virtually unbreakable. For use with the DuPont Technelite® Generator. Magnetic cap keeps the lid on even when inverted. Can be hand held or placed in the single product dispensing station (#600209)

- 600206** The EluTer Elution Vial Shield for DuPont Technelite® Generator (⁹⁹Tc药液瓶后盖式存取罐)
 Ø50×91.5mm N.W:1.82Kg
- 600271** EluTer Top (⁹⁹Tc药液瓶后盖式存取罐上盖)
 Ø38×8mm N.W:0.12Kg
- 600272** EluTer Bottom (⁹⁹Tc药液瓶后盖式存取罐下盖)
 Ø54×15mm N.W:0.33Kg



III.

1. SECURE+™ Unit Dose Syringe Pig

The next generation of the SECURE+™ Unit Dose Syringe Pig has arrived! These tungsten unit dose syringe pigs accept either 3cc or 5cc syringes and weigh only 2.56 lbs. Each pig accommodates a plastic inner sleeve (sold separately) to protect the dose from blood contamination, making each syringe pig a point-of-use sharps container. After use, a red biohazard cap snaps on enclosing the used syringe. Fitted with anti-roll ridges and base for added stability.



- 600268** SECURE+™ Insert Pig Sponge (2.08cm)(100/bag) (⁹⁹Tc药物针管运输罐用黄海绵片)
 Ø20.5mm N.W:0.02Kg/100pcs

- 600297** SECUREplus™ Unit Dose Syringe Shield Pig Foam Insert 3/5cc-12Holes(pair) (⁹⁹Tc药物针管运输罐塑料箱用12孔海绵)
 Top:290×110×25mm N.W:0.18Kg

2. SECURE+™ Unit Dose Syringe Pig Tray

Organize that clutter! The SECURE+™ Unit Dose Syringe Pig Tray allows you to organize your day by holding 12 pigs (#600241) and lids per tray.

- 600241** SECURE+™ Unit Dose Syringe Pig 3/5cc (600218)
 (⁹⁹Tc3cc/5cc药物针管运输罐)
 Ø41.5×199mm Tungsten
 N.W: 1.08Kg
- 600242** SECURE+™ Safety Insert, Tube Body(2310/Case) (⁹⁹Tc针管用塑料管)
 Ø28×122mm 0.08kg/10pcs Plastic
- 600243** SECURE+™ Safety Insert, Long Bio Cap(2720/Case) (⁹⁹Tc针管塑料管长帽)
 Ø30×63mm 0.12kg/10pcs Plastic
- 600244** SECURE+™ Safety Insert, Short Bio Cap(6100/Case) (⁹⁹Tc针管塑料管短帽)
 Ø30×27.5mm 0.10kg/10pcs Plastic
- 600247** SECURE+™ Unit Dose Syringe Pig Tray (⁹⁹Tc药物针管运输罐塑料托盘)
 340×210×50mm Plastic
 N.W:0.55Kg
- 600262** SECURE+™ Unit Dose Syringe Pig, Black Top(10/Pack) (⁹⁹Tc药物针管运输罐塑料上把手)
 Ø35.5×32mm N.W:0.02Kg Plastic
- 600263** SECURE+™ Unit Dose Syringe Pig, Black Bottom(10/Pack) (⁹⁹Tc药物针管运输罐塑料下把手)
 Ø35.5×40mm N.W:0.03Kg Plastic

3. SECURE+™ Unit Dose Syringe Pig 10cc

600261 SECURE+™ Unit Dose Syringe Pig 10cc*Long Bio Cap
(⁹⁹Tc10cc药物针管运输罐)

Ø48.5×223mm Tungsten N.W:1.9Kg

600247 SECURE+™ Unit Dose Syringe Pig Tray
(⁹⁹Tc药物针管运输罐塑料托盘)

Ø340×210×50mm Plastic N.W:0.55Kg

600268 SECURE+™ Insert Pig Sponge(2.08cm)
(100/bag) (⁹⁹Tc药物针管运输罐内黄海绵片)

Ø28.5mm N.W:0.02Kg/100pcs

600275 SECURE+Unit Dose Syringe Shield Pig
Ammo Can Foam Insert 10cc
(⁹⁹Tc药物针管运输罐二孔海绵)

Top:290×110×25mm Bottom:285×100×50mm N.W:0.22Kg/pair



4. SECURE+™ Unit Dose Syringe Pig 10cc

600147 FATCat SECURE plus Safety Insert Pig
(⁹⁹Tc3cc/5cc高剂量药物针管运输罐)

Ø42×199.5mm Tungsten N.W:2.33Kg

600150 FATCat SECURE plus Safety Insert Pig Ammo
Can Foam Insert (⁹⁹Tc3cc/5cc高剂量药物
针管运输罐铁箱用海绵)

Top:280×140×60mm Bottom:280×140×60mm

N.W:0.06Kg

600247 SECURE+™ Unit Dose Syringe Pig Tray
(⁹⁹Tc药物针管运输罐塑料托盘)

340×210×50mm Plastic N.W:0.55Kg



VI. SECURE® Bio Cap Dispenser

You can never find the right Bio Cap when you need it. Now with the New SECURE Bio Cap Dispenser it mounts on the wall right where you need it. It's made of a durable plastic and can hold up to long and /or short Bio Caps in one unit. No more need to look around, they are right on the wall.

600248 SECURE Bio Cap Dispenser

(⁹⁹Tc针管塑料管上帽挂箱)

Dimensions: 502×178×280mm

Material: Plastic

N.W:1.30Kg



V. Hamilton Vial Heating Shield

Made from Tungsten, this Hamilton Vial Heating Shield holds vials the size of 2.65cm in width. Fits within standard aluminum heat blocks. Easy button top and vial cut to lift vials with the use of forceps.

600251 Hamilton Vial Heating Shield

(⁹⁹Tc3cc/5cc药剂瓶蒸馏防护罐)

Size: Ø35×72.7mm

N.W:0.49Kg



VI. Defender™ Dose Drawing Syringe Shields

- *2 mm solid tungsten flange helps shield the hand when withdrawing liquid from a vial.
- *Flange is easily removed to allow transition from drawing dose to patient injection.
- *9 mm thick glass-5.2g/cc gives the greatest protection of any glass in any syringe shield and is easily replaced.
- *Twist Grip-insert any BD® syringe, twist-turn and the syringe is held firmly. No screws, pins or springs to strip out or break-just a simple twist.
- *Easy viewing window -large, optically polished glass allows maximum light entry, and reflective white barrel lining makes viewing syringe contents easy.



- 600287 Latchkey 5cc- Replacement Lead Glass (⁹⁹Tc5cc可视窗铅玻璃)
51.6×11.9×10.4mm N.W.:.03Kg
- 600288 Latchkey 10cc- Replacement Lead Glass (⁹⁹Tc10cc可视窗铅玻璃)
65.6×11.9×9.7mm N.W.:.03Kg
- 600231 Safe Tview™ Syringe Shield 3cc (⁹⁹Tc3cc可视窗尖头式抽药针管防护套)
Ø19.5×63mm N.W.:.0.18Kg
- 600232 Safe Tview™ Syringe Shield 5cc (⁹⁹Tc5cc可视窗尖头式抽药针管防护套)
Ø23×61mm N.W.:.0.22Kg
- 600233 Safe Tview™ Syringe Shield 10cc (⁹⁹Tc10cc可视窗尖头式抽药针管防护套)
Ø25.7×75mm N.W.:.0.31Kg

- 600202 Defender™ Dose Drawing Syringe Shield 3cc* (⁹⁹Tc3cc可视窗园板式抽药针管防护套)
Ø70×63mm N.W.:.0.32Kg
- 600203 Defender™ Dose Drawing Syringe Shield 5cc* (⁹⁹Tc5cc可视窗园板式抽药针管防护套)
Ø70×61mm N.W.:.0.38Kg
- 600204 Defender™ Dose Drawing Syringe Shield 10cc* (⁹⁹Tc10cc可视窗园板式抽药针管防护套)
Ø70×75mm N.W.:.0.43Kg
- 600234 Safe Tview™ 3cc- Replacement Lead Glass (⁹⁹Tc3cc可视窗铅玻璃)
53.7×11.9×9.4mm N.W.:.0.02Kg
- 600235 Safe Tview™ 5cc- Replacement Lead Glass (⁹⁹Tc5cc可视窗铅玻璃)
51.6×11.9×10.4mm N.W.:.0.03Kg
- 600236 Safe Tview™ 10cc- Replacement Lead Glass (⁹⁹Tc10cc可视窗铅玻璃)
65.6×11.9×9.7mm N.W.:.0.03Kg
- 600228 Defender™ 3cc- Replacement Lead Glass (⁹⁹Tc3cc可视窗铅玻璃)
53.7×11.9×9.4mm N.W.:.0.02Kg
- 600229 Defender™ 5cc- Replacement Lead Glass (⁹⁹Tc3cc可视窗铅玻璃)
51.6×11.9×10.4mm N.W.:.0.03Kg
- 600230 Defender™ 10cc- Replacement Lead Glass (⁹⁹Tc10cc可视窗铅玻璃)
65.6×11.9×9.7mm N.W.:.0.03Kg
- 600276 Latchkey Syringe Shield™ 3cc-W/Lead Glass Window (⁹⁹Tc3cc可视窗、勾、椭圆式抽药针管防护套)
Ø19.5×64mm N.W.:.0.19Kg
- 600277 Latchkey Syringe Shield™ 5cc-W/Lead Glass Window (⁹⁹Tc5cc可视窗、勾、椭圆式抽药针管防护套)
Ø23×62mm N.W.:.0.25Kg
- 600278 Latchkey Syringe Shield™ 10cc-W/Lead Glass Window (⁹⁹Tc10cc可视窗、勾、椭圆式抽药针管防护套)
Ø25.7×68.2mm N.W.:.0.33Kg
- 600286 Latchkey 3cc- Replacement Lead Glass (⁹⁹Tc3cc可视窗铅玻璃)
53.7×11.9×9.4mm N.W.:.0.02Kg

VII.

1, Lifetime™ Syringe Shields

*The best shielding syringe shield for patient injection! No lead glass window to crack or fog, just solid tungsten.

*Fully exposed needle hub allows you to visually check for correct venous insertion prior to injection. Permits the use of three-way stopcocks.

*Twist Grip-insert any BD® syringe, twist-turn and the syringe is held firmly. No screws, pins or springs to strip out or break-just a simple twist.

*Uniquely serial numbered solid tungsten syringe shields can be assigned to staff members.

600212	Lifetime™ Syringe Shield 3cc* (⁹⁹ Tc3ccBD罗纹式抽药针管防护套)
	Ø 17×57mm N.W: .0.13Kg
600213	Lifetime™ Syringe Shield 5cc* (⁹⁹ Tc5ccBD罗纹式抽药针管防护套)
	Ø 20.5×52mm N.W: .0.15Kg
600214	Lifetime™ Syringe Shield 10cc* (⁹⁹ Tc10ccBD罗纹式抽药针管防护套)
	Ø 23×70 mm N.W: .0.23Kg
600215	Lifetime™ (Monoject) Syringe Shield 3cc* (⁹⁹ Tc3cc马纳甲罗纹式抽药针管防护套)
	Ø 15.8×61mm N.W: .0.15Kg
600216	Lifetime™ (Monoject) Syringe Shield 5cc* (⁹⁹ Tc5cc马纳甲罗纹式抽药针管防护套)
	Ø 21.2×55mm N.W: .0.17Kg
600217	Lifetime™ (Monoject) Syringe Shield 10cc* (⁹⁹ Tc10cc马纳甲罗纹式抽药针管防护套)
	Ø 24.5×70mm N.W: .0.23Kg

2, Latchkey Syringe Shield

*So simple, it's revolutionary! The only full-length syringe shield that covers the entire length of the syringe barrel and attaches securely with a simple latchkey. Great for unit-dose users!

*Outstanding design. No screws, pins or springs to break or strip out -just a simple twist into the latchkey holds every syringe firmly.

*Fully exposed needle hub allows you to visually check for correct venous insertion prior to injection. Permits the Use of three-way stopcocks.

*Available with or without a window in 3, 5 and 10cc sizes to fit BD syringe. Terumo Syringe Is available in a 2.5cc size.

*Models with windows have a large, optically polished window that allows maximum light entry, and reflective white barrel lining makes viewing syringe contents easy.

--9 mm thick glass-5.2g/cc gives optimum protection and is easily replaced.

--Window is recessed into offset solid tungsten barrel to prevent shine-through.

*Models without windows have no lead glass window to crack or fog, just solid tungsten for maximum shielding -great for unit dose users.



600237	Latchkey Syringe Shield™ 3cc(⁹⁹ Tc3cc尖、勾、椭圆式抽药针管防护套)
	Ø17×68mm N.W: .0.16Kg
600238	Latchkey Syringe Shield™ 5cc(⁹⁹ Tc5cc尖、勾、椭圆式抽药针管防护套)
	Ø21×63mm N.W: .0.23Kg
600239	Latchkey Syringe Shield™ 10cc(⁹⁹ Tc10cc尖、勾、椭圆式抽药针管防护套)
	Ø 23×81mm N.W: .0.26Kg

VIII. SECURE Syringe Shield

*The only syringe shield design that effectively reduces the risk of radiation exposure and blood contamination.

*Exclusive clamshell design opens to allow uncapped used syringes to fall into a waste container instead of backing the syringe back through the syringe shield.

*Overlapping tungsten doors provide optimal shielding.

*Easy viewing window -large, optically polished glass allows maximum light entry, and reflective white barrel lining makes viewing syringe contents easy.

- 600200 SECURE Syringe Shield @3cc
(⁹⁹Tc3cc带可视窗自卸式抽药针管防护套)
N.W: .0.23Kg
- 600279 SECURE Syringe Shield @5cc(ea)
(⁹⁹Tc5cc带可视窗自卸式抽药针管防护套)
N.W: .0.23Kg
- 600201 SECURE- Replacement Lead Glass
(⁹⁹Tc1cc/3cc/5cc针管防护套可视窗铅玻璃)
44×12.2×6.3mm N.W: .0.017Kg



IX.

- 600292
Latchkey Beta/Gamma
Syringe Shield 10cc
(⁹⁹Tc10cc有机玻璃、钨合金、复合
式抽药针管防护套)

Size: Ø31×76mm
Material: Tungsten ,Acrylic
N.W: .0.41Kg



X.

- 600152
Centrifuge Tube Shield Cup
(⁹⁹Tc药剂瓶离心用防护杯)

Size: Ø50×85mm
Material: Tungsten
N.W: .0.78Kg



XI.Tech L-Block

The Tech L-Block is specially designed to protect staff members working with ^{99m}Tc isotopes. This lead-containing L-Block gives full protection to the techologist's torso while the leaded glass window allows the technologist full view of the work area and offers shielding from radiation. This is an excellent product for facilities that compound significant quantities of ^{99m}Tc based radiopharmaceuticals.

600259 Tech L-Block (^{99}Tc 药物分装L型防护屏)

Dimensions: 350×411×587mm

Material: Lead ,Stainless steel

Lead: thick 12mm

Lead Glass: 395×300×28mm

Density of Glass: 5.2cm³

View through the Glass: High quality optical grade

Finish: Smooth stainless steel

Weight: 72 kg



XII. 680004 QC Staition (^{99}Tc 质控用不锈钢支架)

Size: 635×559×298mm

Material :Stainless Steel

N.W:8.5Kg



Xiii. 600140 Qc L-Block (^{99}Tc 低剂量药物分装防护L型防护屏)

Size: 255×255×397mm

Material: Lead, Stainless Steel

Lead: thick 6mm

Lead Glass:

238×208×15mm

N.W: 17Kg



XiV.Tech Clean Air Hood

This drawing station, made of easy-to-clean stainless steel, features a built-in well for the dose calibrator, custom-sized shielded wells for waste collection, a fully shielded base and an integral L-Block style shield. In addition, the drawing station features a 0.3 micron filter, built-in lighting and a pulley and cable assembly system to raise and lower syringes from the ion chamber. This device is the perfect solution for those facilities that use large amounts of activity and prepare doses for many studies.

SPECIFICATIONS:

Working Area: 1100×710×740 Lead Shielded
 Color: Stainless Steel
 Finish: Smooth
 Electrical: 115V, 60Hz, single phase, 15 amps, 9' power cord, 4 wire
 Final Filter: 18.00×48.00×3.00 Aluminum
 Frame, 99.909% efficient, .3 microns w/neoprene air seal gasket w/spring tension, spring tension compression clamps.
 Total CFM: 600
 Lamps: 3-F48T12CW fluorescent cool white, slimline
 Ballast: F36/48T Slimline, 115V, 60Hz, single phase
 Motor: 0.5"hp, 1200 rpm
 Control: Solid state electronic motor speed control
 Material: 16 gauge cold rolled steel
 Work Top: 0.5" formed stainless steel front/rear return panels
 Outlet: 115V, 60Hz, single phase, 5 amps, grounded, universal jack
 Hoist: Pulley and cable assembly for inserting and removing product from ion chamber
 Waste Chute: Lead shielded 8" I.D×12" deep" cylinder
 Shield: Stainless steel personnel shield lead lined bottom and picture frame with 1" lead glass
 Sleeve: Stainless steel formed calibrator holder
 Slide: Stainless steel calibrator slide rail, adjustable for viewing

TABLE STAND

Size: 1280×840×980mm
 Color: Stainless Steel
 Finish: Smooth
 Material: 11 gauge, 2×2 cold rolled steel tube and 16 gauge 1×1 cold rolled steel
 Stand with wheels: movable or stand

Size: 1280×840×2540mm
 Lead : thick 6mm
 Lead Glass (铅玻璃): 510×380×20mm
 Material: Stainless Steel
 Weight: 554kg



600266 Tech Clean Air Hood
 (⁹⁹Tc药物分装100级无菌净化通风柜)

XV.

600253 Tech Lead Waste Barrel (1.27cm)

Tech Lead Waste Barrel(1.27cm)

Surround in Stainless Steel!!!

SPECIFICATIONS:

Height from Ground: 800mm

Handles: 4 every 90 degrees

Wheels: Bolted to removable circular bottom

Wheels Base: 88.9. This lifts barrels 5 inches off the ground.

Alternative wheels are 127 in diam. Height at

Cover handles: 832mm. Height at lid handles: 857mm

600253 Tech Lead Waste Barrel (1.27cm)

(⁹⁹Tc高剂量废弃物存储桶)

Ø466×845mm Material: Lead, Stainless Steel

Lead: thick 12.5mm N.W: 301Kg

600260 Tech Lead Waste Barrel (0.75cm)

(⁹⁹Tc低剂量废弃物存储桶)

Ø466×845mm Material: Lead, Stainless Steel

Lead: thick 7.5mm N.W: 211Kg

600273 Tech Lead Waste Barrel (.5") Top (⁹⁹Tc低剂量废弃物存储桶上盖)

Ø466×65mm Material: Lead, Stainless Steel

Lead: thick 60mm N.W: 45Kg



600264

Dolly: Iodine/Tech Lead Barrels

(⁹⁹Tc废弃物存储桶底盘车)

Ø480×190mm Material: Stainless Steel

N.W : 8Kg

680001

Lead Barrel Lid winch (⁹⁹Tc废弃物存储桶上盖起吊架)

700×750×1248mm Material: Stainless Steel

Square tube(方管): 2×50×50

N.W: 20.5Kg

680002

Red Waste Bucket (⁹⁹Tc废弃物存储桶内红色塑料桶)

Ø380×580mm

Material: plastic

N.W: 1.7Kg

XVI.

600265 Generator Box

This shielded stainless steel case is designed to store multiple generators. The Generator Box completely surrounds the generators, offering 360 degree shielding. It is the ideal solution at those facilities that have an in-house nuclear pharmacy responsible for compounding radiopharmaceuticals.

Size : 1150×600×1055mm

Material: Lead, Stainless Steel

Lead: thick (Side: 9mm)

(top/bottom : 6mm)

Weight: 383kg



⁹⁹Tc母液罐存储箱

Xvii. Tech Clean Air Hood(Small Type)

This drawing station ,made of easy-to -clean stainless steel, features a built-in-well for the dose calibrator, custom-sized shielded wells for waste collection , a fully shielded base and an integral L-Block style shield. In addition, the drawing station features a 0.3 micron filter, built-in lighting and a pulley and cable assembly system to raise and lower syringes from the ion chamber. This device is the perfect solution for those facilities that use large amounts of activity and prepare doses for many studies.



- 700266 Tech Clean Air Hood(Small Type)**
 (⁹⁹Tc药物分装小型100级无菌净化通风柜)
 Size: 120×800×2270mm
 Material: Stainless Steel and Lead
 N.W: 370kg

- 700260 Lead Waste Barrel(Small Type)**
 (⁹⁹Tc低剂量废弃物小型储存桶)
 Size: Ø360×625mm
 Material: Stainless Steel and Lead
 N.W: 123.2kg

- 680005 Casters&Wheels**
 (万向式可升降脚轮)
 Size: 120×110mm
 N.W: 1.42Kg



XVIII.The EluTer Elution Vial Shield for DuPont Technelite Generator (type 2)

These vials shield are made of solid tungsten.



- 600206-2 The EluTer Elution Vial Shield for DuPont Technelite Generator (type 2)**
 (⁹⁹Tc药液瓶后盖式存取罐)
 Ø 44×83.4mm
 N.W: 1.48Kg

- 600206-2 EluTer-2 Body**
 (⁹⁹Tc药液瓶后盖式存取罐)
600271-2 EluTer-2 Top
 (⁹⁹Tc药液瓶后盖式存取罐上盖)
 Ø43.8×8.2mm N.W: 0.15Kg
600272-2 EluTer-2 Bottom
 (⁹⁹Tc药液瓶后盖式存取罐下盖)
 Ø46×11.9mm N.W: 0.22Kg

XIX.The EluTer Elution Vial Shield for DuPont Technelite Generator (type 3)

These vials shield are made of solid tungsten.



600206-3 The EluTer Elution Vial Shield for DuPont Technelite Generator(type R)
 (R-⁹⁹Tc药液瓶后盖式存取罐)
 Ø44×83.4mm
 N.W:1.48Kg

XX.R--Lead Unit Dose Syringe Pig



The unit dose syringe pig is made of lead.It can accepts either 3cc or 5cc syringe.

700241 R--Lead Unit Dose Syringe Pig
 (R-⁹⁹Tc3CC/5CC药物针管铅运输罐)
 Ø37×215.5mm
 N.W:1.36Kg

XXI.R--Lead Magnetic Vial Shield

This Lead Vial Shield assure safe shielding of vials containing radioactive liquids.


700208 R--Lead Magnetic Vial Shield
 (R-⁹⁹Tc药液长瓶铅存取罐)
 Ø64×91mm
 N.W:2.5Kg





XXII.

- 600231-1cc Safe Tview Syringe Shield 1cc-BD1ml Luer Lock**
 (⁹⁹Tc1cc带窥视窗尖头式抽药针管防护套)
 Ø19.5×71mm
 Material: Tungsten
 N.W:0.23Kg
- 600231-2cc Safe Tview Syringe Shield 2cc-Braun 2ml**
 (⁹⁹Tc2cc带窥视窗尖头式抽药针管防护套)
 Ø19.5×48mm
 Material: Tungsten
 N.W:0.14Kg

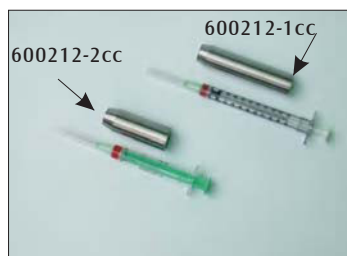

- 600212-1cc Lifetime Syringe Shield 1cc-BD1ml Luer Lock**
 (⁹⁹Tc1cc罗纹式抽药针管防护套)
 Ø16.5×71mm
 Material: Tungsten
 N.W:0.17Kg
- 600212-2cc Lifetime Syringe Shield 2cc-Braun 2ml**
 (⁹⁹Tc2cc罗纹式抽药针管防护套)
 Ø17.5×48mm
 Material: Tungsten
 N.W:0.11Kg



XXIII.

- 700152 R-Lead Centrifuge Tube Shield Cup**
 (R-⁹⁹Tc药剂瓶离心用铅防护杯)
 Ø95×158mm
 Material: lead
 N.W:2.8Kg

Iodine

I. 600108 Iodine Dispensing System

The Iodine Dispensing System permits dispensing of radioiodine without having to hold the pig. The 2.5cm thick solid tungsten container holds a 10ml vial and safely shields up to Curie quantities of iodine ¹³¹I. With its tungsten "cork" magnetically held in place, it is designed to be held inverted in the dispensing station to permit safe compounding of iodine capsules or solutions. The stainless steel base provides shielding when standard lead bricks are in place. Its fully adjustable armature permits positioning to any height or angle.

The complete 600108 Iodine Dispensing System* includes:

- 600109 Iodine Dispensing Pig* (¹³¹I药液瓶存放罐)
Ø77×117mm Material: tungsten N.W:7.4Kg
- 600110 Iodine Syringe Shield 0.5/1cc
(¹³¹I1cc开口、圆板式抽药针管防护套)
Ø70×59mm Material: tungsten N.W:0.44Kg
- 600111 Stainless Steel Base, Armature & Ring
(抽药用支架)
307×210×60mm Material: Stainless Steel N.W:5Kg
- 600132 Nordion 3ml V-Vial Adapter
(¹³¹I药物存放罐内衬杯)
Ø25.5×48mm Material: Plastic N.W:0.02Kg



II.

1. 600100 Piglet™ (¹³¹I吸碘瓶防护罐)

This solid tungsten delivery container provides optimal shielding for iodine I-131 with minimum weight (only 4.46Kg). Its direct threaded lid assures safe containment of either liquid or capsule dosages. The Piglet™ is capable of safely shielding high activity carcinoma ablation doses.

2. 600125 Piglet™ Sling (¹³¹I吸碘瓶防护罐尼龙袋)

The ideal accessory for the Piglet™ or Piglet2™. Make sure they are the best-dressed Piglet™ in town! These color-coded slings are crafted from ballistic nylon and can help keep exposure ALARA while facilitating easy handling.

- 600100 Piglet™ (¹³¹I吸碘瓶防护罐)
Ø60×96.4 Tungsten N.W:4.6Kg
- 600125 Piglet™ Sling (¹³¹I吸碘瓶防护罐尼龙袋)
N.W:0.08Kg/10pcs
- 600102 Piglet™ Ammo Can Foam Insert
(¹³¹I防护罐铁箱海绵)
280×105×105 N.W:0.07Kg
- 600116 Piglet™ Lid Gasket-each (¹³¹I防护罐F4垫圈)
Ø59×Ø40.9×0.8 N.W:0.02Kg/10pcs
- 600107 Charcoal Donut Piglet™-100/bag (¹³¹I防护罐内活性炭吸水垫片)
N.W:0.001Kg/10pcs



III.

1. 600101 Piglet2™ (¹³¹I检测瓶防护罐)

Piglet2™ is the smaller cousin to the Piglet™. It is also a solid tungsten delivery container that provides optimal shielding for iodine I-131 with minimum weight (only 2.75 Kg). Its direct threaded lid assures safe containment of either liquid or capsule dosages. The Piglet 2™ is capable of safely shielding hyperthyroid or whole body scan doses.


2. 600124 Piglet™ Pooper (¹³¹I检测瓶开瓶启子)

The Piglet™ Pooper is a remote-handling tool that helps remove and replace the snap cap lid within the Piglet™ or Piglet 2™. This stainless steel tool is a must for every department and can help keep exposure ALARA. There is no more need to physically touch the vials because of the easy small tension grip used for opening and closing.

3. 600120 Gotti Cup (¹³¹I检测瓶抽药用防护杯)

The Gotti Cup is constructed of tungsten and is designed for use in I-131 NaI capsule compounding. Using the Gotti Cup dramatically increases the level of shielding available during the capsule compounding process.

- 600101 Piglet 2™ (¹³¹I检测瓶防护罐)
Ø 51.8×88.6mm N.W: 2.84Kg
- 600120 Gotti Cup (¹³¹I检测瓶抽药用防护杯)
Ø 33×37.8mm N.W: 0.53Kg
- 600126 Piglet 2™ Sling (¹³¹I检测瓶防护罐尼龙袋)
N.W: 0.07Kg/10pcs
- 600117 Piglet 2™ Lid Gasket-each
(¹³¹I检测瓶防护罐F4垫圈)
Ø 49.5×Ø 40.8×0.8
N.W: 0.01Kg/10pcs
- 600124 Piglet™ Pooper (¹³¹I检测瓶开瓶启子)
N.W: 0.05Kg
- 600123 Charcoal Donut Piglet 2™-100/bag
(¹³¹I防护罐内活性炭吸水垫片)
N.W: 0.001Kg/10pcs

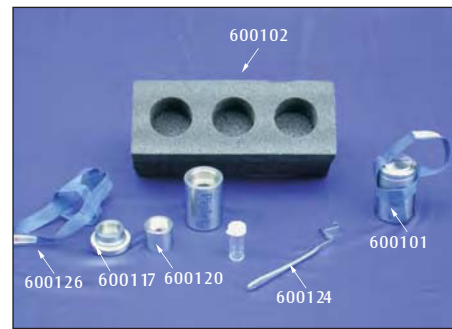

- 600102 Piglet™ Ammo Can Foam Insert (¹³¹I防护罐铁箱海绵)
280×105×105mm
N.W: 0.07Kg


IV.



- 600129 Piglet 3™ (¹³¹I吸碘方瓶防护罐)
Ø 85×119mm Tungsten N.W: 10.7Kg
- 600134 Piglet 3™ Ammo Can Foam Insert
(¹³¹I方瓶防护罐铁箱单孔海绵)
280×140×120mm N.W: 0.13Kg
- 600135 Piglet 3™ SECURE+™ Ammo Can Foam Insert
(¹³¹I方瓶防护罐塑料箱双孔海绵)
270×140×90mm N.W: 0.08Kg
- 600136 Piglet 3™ Sling (¹³¹I方瓶防护罐尼龙袋)
N.W: 0.2Kg/10pcs
- 600145 Piglet 3 Round Circular Sponge
(50/bag)
(¹³¹I方瓶防护罐黄海绵方孔环)
N.W: 0.006Kg/10pcs (D34 孔28*28)
- 600146 Piglet 3 Square Sponge (50/bag)
(¹³¹I方瓶防护罐内黄海绵垫28×28)
N.W: 0.004Kg/10pcs

V.

- 600284 SECUREplus™ Delivery Multi Dose Nylon Case Foam Insert -one Hole
(ϕ 56-一孔塑箱海绵)
105×95×77mm N.W:0.02Kg
- 600285 SECUREplus™ Delivery Multi Dose Nylon Case Foam Insert -Two Hole
(ϕ 56-二孔塑箱海绵)
140×105×77mm N.W:0.02Kg
- 600148 Multi Dose Ammo Can Foam Insert(10/bag)
(ϕ 56-二孔铁箱海绵)
140×75×77mm N.W:0.02Kg


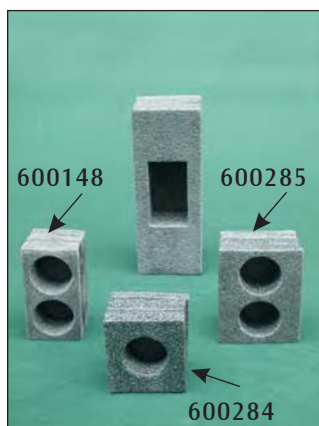


VI.

- 600400 Seed Dipper (^{125}I 粒子检测器)
Size: $\phi 90 \times 290\text{mm}$
Material: Acrylic
N.W: 0.11Kg

VII.

- 600130 Iodine Charcoal Filter Box (Sm.)
(^{131}I -蒸气吸附过滤器--7块过滤板)
552×350×310mm
Material: Lead Stainless Steel
N.W: 26.5Kg
- 313673 Charcoal Air Filer Type 787
(I-131) For Iodine Charcoal Filter Box & Glove Box (^{131}I -蒸气吸附板)

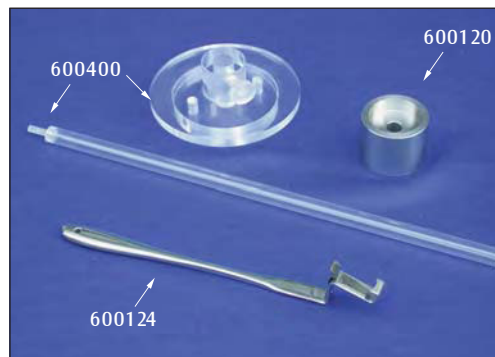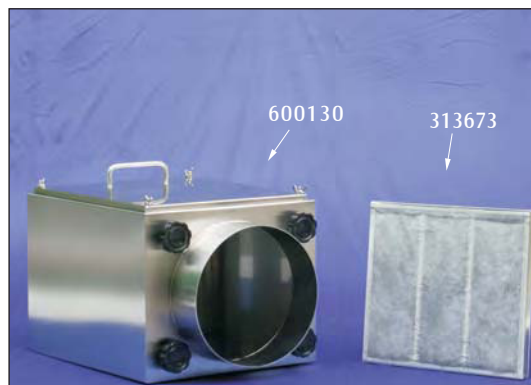

VIII.

- 600155 I-131 Waste Storage Cabinet (^{131}I 废弃物储存箱)



Size: 1312×636×490mm
Material: Lead
Stainless Steel
Lead: thick side 9mm
Top, Bottom 6mm
N.W: 443.5Kg

IX.

600128 Iodine Lead Waste Barrel
Barrel (1.9cm thick)
 (¹³¹I高剂量废弃物存放桶--PET废弃物可用)

Iodine Lead Waste Barrel
 Surrounded in Stainless Steel!!

SPECIFICATIONS:
 Barrel: 470mm diam
 Height from Ground: 800mm
 Handles: 4 every 90 degrees
 Wheels: Bolted to removable circular bottom
 Size: Ø466×845mm
 Lead: thick 19mm
 Material: Lead, Stainless Steel
 N.W: 338Kg



- 600128 Iodine Lead Waste Barrel**
(1.9cm thick)
 (¹³¹I高剂量废弃物存放桶--可用于PET废弃物)
 Ø466×845mm N.W: 338Kg
- 600141 Iodine Lead Waste Barrel**
Top(¹³¹I高剂量废弃物存放桶上盖)
 Ø466×65mm N.W: 55Kg
- 600264 Dolly: Iodine/Tech Lead**
Barrels (¹³¹I/ ⁹⁹Tc废弃物储存桶移动车)
 Ø480×190mm Stainless Steel
 N.W: 8.0Kg
- 680001 Lead Barrel Lid winch**
 (¹³¹I高剂量废弃物存放桶上盖起吊架)
 700×750×1248mm
 Material: Stainless Steel
 N.W: 20.5Kg
- 680002 Red Waste Bucket**
 (¹³¹I/ ⁹⁹Tc废弃物存放桶内红色塑料桶)
 Ø380×580mm
 Material: Plastic
 N.W: 1.7Kg



X.I-131-L-Block (¹³¹I药物分装L型防护屏)

This shield is designed for facilities that prepare and dispense high activity ¹³¹I therapy doses.

The shield offers excellent protection from beta and gamma radiation and has a fully shielded floor. The lead glass is 3.5 inches thick to shield out the high energy gamma photons from this radionuclide.

600133 **I-131-L-Block** (¹³¹I药物分装L型防护屏)

SPECIFICATIONS:

Dimensions: 350×310×448mm

Lead: thick 36mm

Lead Glass (铅玻璃): 295×228×68mm

Density of Glass: 5.2cm³

View through the Glass: High quality optical grade Finish:

Smooth stainless steel

Material: Lead,Stainless Steel

Weight: 102.7kg



XI. 600131 Glove Box & Table (¹³¹I药物分装手套箱)

Size: 1200×800×2175mm

Material: Lead,Stainless steel

Lead (thick): front(前面) 20mm
back(后面) 50mm

24 WKg



XII.NeRD[®]

*This innovative needle recapping device facilitates OSHA compliant one-handed Needle recapping.

*The suction cup, or permanent adhesive base provides stable positioning and it accepts most needle brands sizes and gauges.

*Its small size makes it easily transportable and it Won't clutter the work area.

*Made of solid rubber with a highly visible color.

*It is easily cleaned and decontaminated.



600119 NeRD (带吸盘针帽座)
N.W:0.03Kg

XIII.Tweezers - 9.75 Inch

No more need to physically touch chromatography strips. Surgical grade steel. Great for use in quality control procedures.

600114 Tweezers - 9.75 Inch (不锈钢镊子)
N.W:0.07Kg

These stainless steel forceps are excellent for handling small items in the hot lab. Use the

XIV.Forceps - 10 Inch

inverse square law to your advantage!

600115 Forceps 10",Curved (止血钳-弯有槽)
N.W:0.07Kg

600138 Forceps Guard(round) (¹³¹I操作钳圆形防护屏)
N.W:0.08Kg

600139 Forceps Guard(trapezoid) (¹³¹I操作钳方形防护屏)
N.W:0.1Kg

600162 Forceps 10",Straight (止血钳-直有槽)
N.W:0.06Kg



XV.R--Iodine Pig 1 (R-¹³¹I吸碘瓶防护罐)



This solid tungsten delivery container provides optimal shielding for iodine I-131 with minimum weight (only 5.43Kg). Its direct threaded lid assures safe containment of either liquid or capsule dosages. The Iodine Pig 1 is capable of safely shielding high activity carcinoma ablation doses.

600100-2 R--Iodine Pig 1 (R-¹³¹I吸碘瓶防护罐)
Material: Tungsten
Size: $\varnothing 66 \times 104\text{mm}$
N.W: 5.43Kg

XVI.R--Iodine Pig 2 (R-¹³¹I检测瓶防护罐)



Iodine Pig 2 is the smaller cousin to the Iodine Pig 1. It is also a solid tungsten delivery container that provides optimal shielding for iodine I-131 with minimum weight (only 3.64Kg). Its direct threaded lid assures safe containment of either liquid or capsule dosages. The Iodine Pig 2 is capable of safely shielding hyperthyroid or whole body scan doses.

600101-2 R--Iodine Pig 2 (R-¹³¹I检测瓶防护罐)
Material: Tungsten
Size: $\varnothing 57.7 \times 94\text{mm}$
N.W: 3.64Kg

XVII. 600151 Squeezer QC Syringe Shield (I-¹³¹药物针管W-橡胶防护套)

Material: Rubber
Size: $\varnothing 17 \times 68\text{mm}$
N.W: 0.06Kg



XVIII. 600502 Lead Loose Seed Holder (I-¹²⁵粒子瓶铅防护罐)

Material: Lead
Size: $\varnothing 28 \times 66.6\text{mm}$
N.W: 0.30Kg



Stainless Steel

I.

- 600408 SS Cabinet:(A)-24" Drawers (24"抽斗立橱)
610×610×890mm Material: Stainless Steel
N.W: 60.3Kg
- 600409 SS Cabinet:(C)-36" Cabinet(36"双门双抽立柜)
915×610×890mm Material: Stainless Steel
N.W: 62.5Kg
- 600410 SS Cabinet:(C)-48" Cabinet(48"双门双抽立柜)
1220×610×890mm Material: Stainless Steel
N.W:75.8Kg
- 600411 SS Cabinet:(D)-48"Sink (48"带水斗双门立柜)
1220×610×890mm Material: Stainless Steel
N.W:50.2Kg



II.

- 600412 SS Cabinet: (E) -24"Overhead Open Cabinet (24"有门挂柜)
610×356×787mm Material: Stainless Steel N.W:33.5Kg
- 600413 SS Cabinet: (F) -48"Overhead Open Cabinet (48"有门挂柜)
1220×356×787mm Material: Stainless Steel Glass :582×460×5mm N.W:59Kg



III.

600414 SS Cabinet: (G) -48"Overhead Open Self
(48"无门挂柜)

Size: 1220×356×787mm Material: Stainless Steel N.W:41Kg



IV.

600422 Conveyer Table (不锈钢四方桌)

Size: 990×660×984mm
Material: Stainless Steel
N.W:32Kg



V.

600423 Table Base for 2 Waste Barrels (二桶用桌架)
1372×762×930mm Material: Stainless Steel N.W:27Kg

600423 Table Base for 2 Waste Barrels (Special) (一大一小桶用桌架)
1676×762×930mm Material: Stainless Steel N.W:29Kg



VI. **600424 Table Base for 3 Waste Barrels (三桶用桌架)**

1981×762×930mm Material: Stainless Steel N.W: 33Kg



VII. **600427 Lead Funnel & Tube (放射性废弃物防护套)**

Size (1): 12×Ø160×Ø137×160mm
Size (2): 12×Ø190×105mm
Material: Lead, Stainless Steel
Lead: thick 8mm
N.W: 11.8Kg



VIII. **600416 SS Cabinet:(SL)-Right Support Leg** (不锈钢台板右面支撑脚)

660×170×889mm
Material: Stainless Steel
N.W: 3Kg

600417 SS Cabinet:(SL)-Left Support Leg (不锈钢台板左面支撑脚)

660×170×889mm
Material: Stainless Steel
N.W: 3Kg



IX.

600418 SS Cabinet:(CG)- Corner Guard (不锈钢转角)

63.5×63.5×1830mm
Material: Stainless Steel
N.W: 3Kg

Others

I.

1. 600256 Ashley Dose Calibrator Funnel Dipper

Imagine a dose calibrator dipper where you don't have to thread the needle! The unique funnel design positions the syringe in place, you just toss it in. Accommodates 0.5,1,3,5,6,10,1220 and 60 ml syringes!

SPECIFICATIONS:

Construction: Acrylic

2. 600257 Standard Dose Calibrator Dipper

This dose calibrator dipper is specially designed to hold syringe and vials. Fits in all Capintec dose calibrators.

SPECIFICATIONS:

Construction: Acrylic

3. 600258 Well Liner Sleeve

This well chamber insert fits all models of Capintec dose calibrators. The sleeve protects the dose calibrator from contamination and is easily removed for cleaning.

SPECIFICATIONS:

Construction: Acrylic



- | | |
|--------|--|
| 600256 | Ashley Dose Calibrator Funnel Dipper
(活度计配置的有机玻璃测量支架--漏斗型)
N.W:0.06Kg |
| 600257 | Standard Dose Calibrator Dipper
(活度计配置的有机玻璃测量支架--圆型)
N.W:0.10Kg |
| 600258 | Well Liner Sleeve
(活度计有机玻璃防护筒)
N.W:0.23Kg |

II. Hot Box

The Hot Box is a radiopharmaceutical heating containment system designed to contain airborne radioactive contamination that might result during the routine heating steps required in the preparation of certain technetium ^{99m}Tc labeled radiopharmaceuticals. These radiopharmaceuticals might be, but are not limited to, any of the following:

- ^{99m}Tc Sulfur Colloid, ^{99m}Tc Mertiatide (MAG3[®]),
- ^{99m}Tc Sestamibi (Cardiolite[®] or Miraluma),
- ^{99m}Tc Apcitide (AcuTect[®]) or ^{99m}Tc Depreotide (Neo Tect[®])



The Hot Box is designed to house a dry heater. This heater must be smaller than 17 cm high, 23cm wide and 30 cm deep to fit within the Hot Box. This dry heater should be equipped with one or two blocks (as appropriate) and a thermometer. The Hot Box is not designed for use with liquid heating baths that utilize water or other low boiling temperature liquids.

The Hot Box is not designed to shield the user from radioactive exposure inherent in the handling of radioactive materials. It is constructed from stainless steel and normal window glass, not lead, tungsten or leaded glass.

600240 Hot Box (^{99}Tc 药物加热器防爆盒)
Size: 327×260×200mm
Material (材料): Lead, Stainless Steel
Lead: thick 6mm
N.W: 27.9Kg

III.SS WBC Case

These stainless steel cases are specially designed for the transport of blood products for labeling procedures .The cases feature the following:

- * All stainless steel construction.
- * Sliding panel door for access (not a hinged door) that can be secured with a security tie for delivery.
- * No foam insert used (nothing to get contaminated with blood)
- * Comes labeled with "This end up" and "Biohazard-Patient Blood Sample "labels already affixed.
- * Door is reversed by customer (exposing "Biohazard" label) to indicate that case is to be returned to pharmacy.
- * Holds two 60 cc syringes or centrifuge tubes (in biohazard Ziploc bags).

600283 SS WBC Case (采血箱)



Size: 375×201×115mm
Material: Stainless Steel
Case: thick 1.2mm
Panel door (盖): thick 2mm
N.W: 3.1Kg

IV.SEED

600401 Seed Mick Cartridge Holder (I¹²⁵粒子不锈钢多孔罐)

Size: Ø63.5×83mm
Material: Stainless Steel
N.W: 0.72Kg

600402 Loose Seed Holder (I¹²⁵粒子不锈钢单孔罐)

Size: Ø29×61.6mm
Material: Stainless Steel
N.W: 0.17Kg



V.



680008 Single Core Filter
(单芯过滤器)
Ø97×58mm
N.W:4.2Kg

VI.

- 680009 J--Tungsten Syringe Shield**
(J-针管防护套)
Ø42×80mm
Material:Tungsten
N.W:1.16Kg
- 680010 J--Tungsten Stopper Shield For Cup Shield**
(J-放射性药液防护杯的二节式防护塞)
Top:Ø20×25mm Bottom:Ø66×104mm
Material:Tungsten



VII.



680011 Tungsten Funnel
(钨合金漏斗)
Ø50×64mm
Material:Tungsten
N.W:2Kg

VIII.

- 680006 W-Cu Alloy Hollow Tube**
(铜钨合金开口管)
Size:Ø 16.5×Ø12×69.5mm
Material:W-Cu Alloy
N.W:0.09Kg
- 680007 W-Cu Alloy Solid Bar**
(铜钨合金实心棒)
Size: Ø15×65mm
Material:W-Cu Alloy
N.W:0.2Kg

